




# HUR Offload Unweighing Device

VERSATILE SOLUTIONS FOR REHAB AND CONDITIONING

# HUR Offload

The HUR Offload is an unweighing system powered by air-pressure that provides adjustable support to patients while they participate in walking, as well as a full range of other therapeutic activities. By reducing body weight, patients are able to engage early in the rehabilitation process in a safe and secure environment.

HUR's unique two-rope system provides enhanced support and safety, while still enabling the ability to turn sideways or backwards. The open design of the HUR Offload provides therapists full access to the patient's lower extremities, and can be used with any treadmill with a maximum width of 950mm. For best results, the HUR treadmill is recommended.

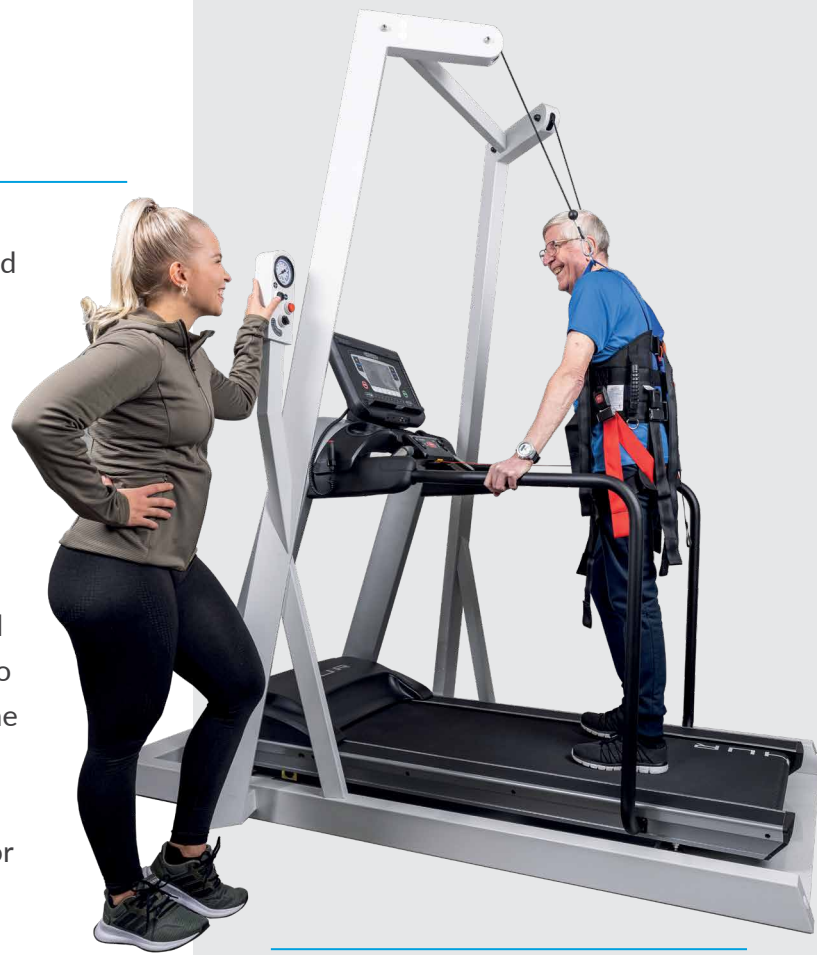

## FEATURES

- A. Air-Pressured Technology**  
Unweighing is achieved through using air pressure (pneumatic) providing constant support in a smooth movement
- B. Affordable Durable Design**  
A high quality and cost-effective solution made in Finland
- C. Incremental Adjustments**  
Small incremental adjustments allows for controlled progress
- D. Easy Access**  
An open design removes barriers and is ideal for older adults and wheelchair users
- E. User Friendly**  
Easily adjust the compensation weight with a simple push of a button
- F. Two-Rope Support System**  
Enhanced support for walking efficiency and stability

### Dimensions:

Max User Weight: 250kg (550 lb)  
Max Unweighing Capacity: 139kg (306 lb)  
Max Width: 950mm (37 in)

### Available in two different heights:

Product 6100 -  
142 x 216 x 239 cm (55.8" x 84.9" x 93.9")  
Product 6200 -  
142 x 216 x 258 cm (55.8" x 84.9" x 101.7")

### WHAT IT IS USED FOR:

- Stroke
- Spinal Cord Injuries
- Vestibular Disorders
- Neurological Disorders
- Balance Training
- Amputees
- Orthopedic

## HUR Unweighing Harness

The HUR Unweighing Harness provides maximum support with optimal comfort. The design equally distributes pressure through the hip straps, gluteal straps and the lumbar-thoracic region without limiting the range of motion in the lower extremities.



## HUR Rehab Treadmill

The HUR Rehab Treadmill is a versatile solution for rehabilitation, physical therapy, and sports conditioning. A zero starting speed, extended parallel handrails, and a low step up height (14cm / 5.5in) makes the HUR Treadmill ideal for a wide range of abilities. In addition, a reversing belt and orthopedic suspension system make it a great tool for walking and other therapeutic activities.



## FEATURES

- A. Zero starting speed**  
A starting speed of zero makes the HUR Rehab Treadmill ideal for rehabilitation environments
- B. Easy Access**  
Low step up height (14cm / 5.5") and extended parallel handrails make getting on and off the treadmill easier
- C. Small, Incremental Adjustments**  
A speed range of .16 - 19 km/h (.1-12 mph) and adjusts in .1 increments
- D. VFX Shock Absorption System**  
An industry leading shock absorption system that is five-times softer than grass
- E. Therapist Features**  
Hand-held remote stop-start switch allows therapists to remain in control of treadmill operation
- F. Rehabilitation Control Panel**  
Multiple heart-rate control programs, fitness tests, and motivational display screens
- G. Reversing Belt**  
Ideal for downhill and retrowalking therapies

### Dimensions:

Max User Weight: 180kg (400lbs)  
Product 6700 - 50 x 147cm (20" x 58")





# YOUR PARTNER FOR **LIFELONG STRENGTH**

---

HUR is dedicated to contributing where it is best = providing the best exercise solutions for preventive, rehabilitative and maintenance training - so that you and your customers can keep on performing miracles every day.



Ab HUR Oy Patamäentie 4, 67100 Kokkola, Finland/Tel +358 6 832 5500 / sales@hur.fi

[www.hur.fi](http://www.hur.fi)